

Fiber Optic ATR-Probes for Harsh Environment



FlexiSpec[®]

High throughput in selected parts
of Near & Mid InfraRed-spectra

ATR-tip shaped for immersion
in liquid flow without dead zone

Robust for industrial applications
in harsh environment up to 250°C

Resistant to high pressure

Compatible with all spectrometers
and automated process-interfaces

FlexiSpec[®] probes from **art photonics** include
new type of Attenuated Total Reflection (ATR-) probes
designed for harsh application conditions.

HT-ATR immersion IR-fiber probes can be used for
process-spectroscopy in Near & Mid IR range to monitor
reactions in-line in a broad temperature range from -100° to
+250°C. They can resist to high pressures up to 200Bar and
used with FTIR or any other IR-spectrometers and spectral
sensors in automated process control with process-interfaces

Applications:

- Remote Reaction Monitoring in-line
in temperature range -100/+250°C
- PAT applications in lab, pilot plant or
industry with process-interfaces for
automated process control
- Polymerization Process Control
- In-situ IR-Spectroscopy for PAT in
Chemical, Petrochemical, Atomic,
Biopharmaceutical & Food Industry




원우시스템즈
WONWOO SYSTEMS

Tel (02) 533-6720, Fax (02) 3289-1293

서울시 동작구 신대방1가길 38 (신대방 719 동작상떼빌) 106동 209호

www.wonwoosystem.co.kr

Transmission spectra of ATR probes. Probe length 1.5m

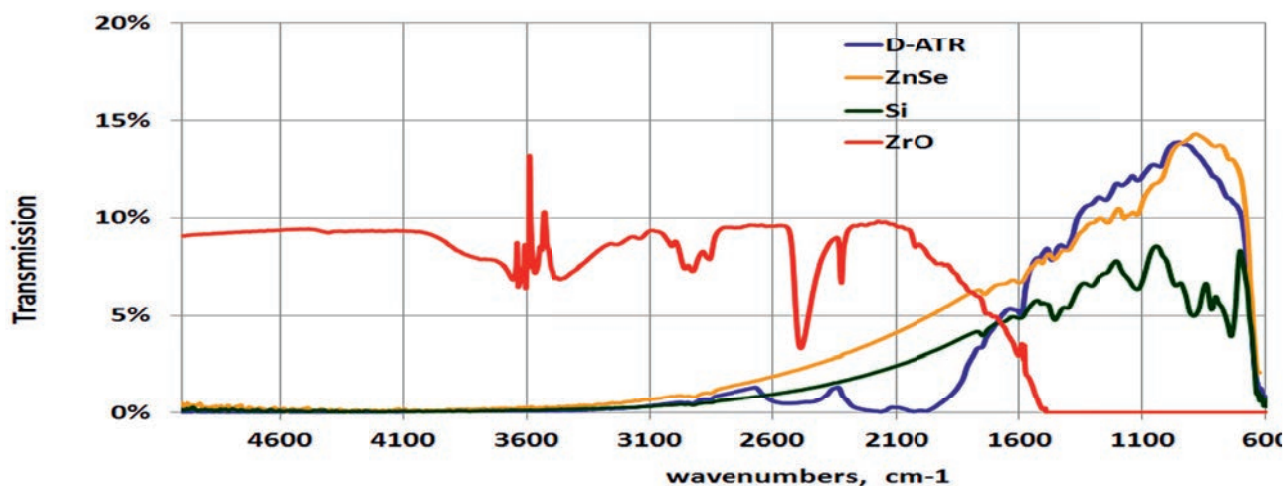

Specification of IR-Fiber HT-ATR-Probes for High Temperatures *FlexiSpec*®

Probe type	Diamond ATR	Silicon ATR	Cubic Zirconium ATR
Transmission range	5.2-17μm (600-1900cm ⁻¹) + 3.2-4.5μm (2300-3100cm ⁻¹)	5.2-17μm (600-3100cm ⁻¹)	1.5-6.5μm (1550-6650cm ⁻¹)
Fiber type	PIR-900/1000 Silver Halide	PIR-900/1000 Silver Halide	Chalcogenide glass (As-S)
Temperature range	-150°C / + 250°C	-150°C / + 250°C	-150°C / + 200°C
Pressure (max)	200Bar	100Bar	100Bar

Common Parameters of HT-ATR-Probes for High Temperatures *FlexiSpec*®

Total Length	1.5m (opt.: 1m to 5m)*
Shaft Length	300 mm (opt.: 300-700 mm)*
Shaft Diameter	12mm
Shaft Material	Hastelloy C22
Length of Legs	500mm (opt.: 200mm to 500mm)
Protective Tube Material	Liquid Tight SS-Conduit, KOPEX-Tube
Minimal Bending Radius	130mm
Input / Output Connectors	Long SMA (opt.: any other type)
Cooling Air Flow Parameters	Excess pressure 0.5Bar, Flow 2300l/h
Inner temperature control	Inside the shaft to control ATR tip temperature

